


Quick Start Guide for the Dyna-Chi Ionic Body Balancer™

1. Read this Quick Start Guide before using the machine. Complete User's Manual available online at www.maxdetox.com.
2. Use a **PLASTIC** tub approximately 17" x 12" x 6" (or higher) from Wal-Mart, Target, etc. The Sterilite and Rubbermaid brands work well, but any plastic tub of this size will suffice. Liners are optional, but the toxic material *will stain the tub*. Professionals use liners, which in the industry are nothing more than trash bags. **This technology does NOT require any special tub**. Metal tubs (including copper tubs) will destroy your power supply.
3. Stand the the array in the middle of the tub. **Do NOT pick the array up by the cables**, this can cause the terminals to loosen. (*Owner's Manual pg. 10-11*)

CORRECT PLACEMENT

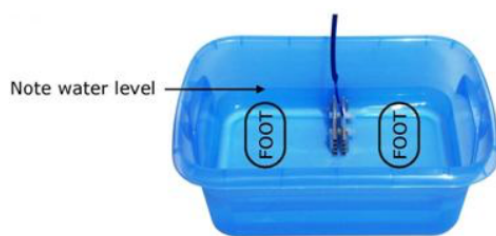

WRONG PLACEMENT



4. Place 1 teaspoon of salt in the tub and fill with **hot** water (as hot as can be tolerated) to the top of the array's bars (*see picture above*). Within minutes, vigorous fizzing should take place around the array. If fizzing does not occur, add additional salt (around the array) by ½ teaspoon increments, not to exceed 3 teaspoons total. **If the power supply becomes hot to the touch you have used too much salt. and should change water.** Note: *This technology does not work well using cold water.*
5. **Plug the power supply into a surge protector** and place the unit **AWAY** from the tub so it cannot fall into the water. Dropping appliances connected to the wall into water can be fatal. The current coming from wall is called "Alternating Current" or "AC". The current generated by the power supply is "Direct Current" or "DC", which is harmless.
6. Plug the gold cable connectors into the color coded ports on the power supply. **Do not remove the plugs by pulling on the cables. Do not submerge the cables.**
7. Place the tub containing the water, salt, and array so that your feet may rest in it comfortably while you are seated. You will also want to have an old towel on hand to dry your feet after your session. The toxins released from cleanse will stain fabric and carpet.
8. Turn on the power supply, sit back and relax, enjoying your treatment.
9. Reminder: If the power supply becomes hot to the touch, you have too much salt in the water! If you feel a tingling around the water line, you may have used too much salt. You are **not** being shocked or harmed if this sensation occurs. Try adding one cup hot water to dilute the tub water and relieve any discomfort, or change water.
10. After use, clean your array PROPERLY. Directions for cleaning may be found at www.ionicbalancer.com/clean. Failure to clean properly will substantially shorten the life of the array.

NOTE: This machine causes a substantial amount of acid to be released from the body. In a small percentages of users, this may cause a temporary rash on feet and ankles following a treatment. If this occurs, discontinue use, wash feet with soap and treat with calamine lotion until the condition clears.

Color of Water	Area of Body	Symptoms
Yellow-Green	Kidneys, Bladder, Urinary Tract, Reproductive Organs	Gynecological, or Prostrate Disorder, Kidney or Bladder Infections
Black Flecks	Heavy Metals	Memory Problems, Chronic Fatigue, Mental "Fog"
Orange	Joints	Arthritis, Rheumatism, Gout, Neuropathy, Fibromyalgia
Brown	Lungs, Cellular Debris	Smoking, Exposure to Second-Hand Smoke, or other forms of Air Pollution
White (Cheese-like Particles)	GI Tract, Any area effected by Yeast Infections	Candida, Migraine, Insomnia, Headaches, Yeast Infections
Dark Green	Gall Bladder, Liver, Stomach	Gall Stones, Weak or Compromised Immune System, Acid Reflux, Nausea
Red Flecks	Blood Clot Material	Arterial Blockage, Bruising, Hypertension, Cold Hands/Feet
White Foam	Lymphatic System	Menses Disorders, Lymph Nodes Tender or Clogged, Skin Disorders or Allergies
Black	Liver	History of heavy drinking or drug use (legal or illegal), Elevated liver enzymes on bloodwork

Warnings:

Do NOT bump or hit the red and black cable terminals on the front of the power supply unit. These can easily break off, which is NOT covered under the warranty.

Use a Surge Protector with the machine to protect it from power surges. Analog machines can easily be damaged by even minor power surges. **Not using a surge protector VOIDS the warranty.**

NOTE to users in North America: With analog power supplies, the higher the saline content the more power is pulled from the power supply.

This technology is affected by body mass. A small person cannot comfortably tolerate as much power as a large person, and therefore needs to add **LESS** salt to the water. **If you feel a biting, uncomfortable tingling when using the machine, you are likely using too much salt in the water. Too much salt added to the water can burn out your power supply which VOIDS the warranty.**

Do NOT use with a metal tub, doing so will cause the power supply to short circuit and burn it out, which is NOT covered under your warranty.

Machine is not recommended for use by children under 100 pounds, or under the age of 16.